


## Performance You Can Feel Results You Can Trust



### Lightweight

The ergonomic design of the lightweight Rainin multichannel Pipet-Lite XLS+ delivers extraordinary comfort and performance, even after hours of pipetting. The feel is balanced and substantial with a fit and finish that is truly world class.

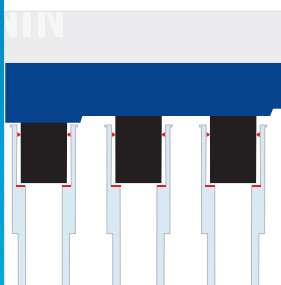

### Lightest touch

Nothing beats the Rainin LiteTouch System for ergonomics and ease of use. The patented LTS shaft and tip design with small seal and positive stop provides a fast, perfect seal across all channels. With LTS and progressive ejection, tips eject easily and with less force.



### Greater control

With an unparalleled sense of control and feedback, Pipet-Lite XLS+ is the easiest and most responsive multichannel pipette you'll ever use. Plunger stroke is smooth and absolutely consistent from start to finish.



### Superior performance

Consistent sample loading across all channels, secure seals, easy tip mounting and the lowest tip-ejection forces are all hallmarks of Rainin multichannel pipettes.



### Pipet-Lite XLS+ Multichannels Unmatched Precision and Performance

The Rainin multichannel Pipet-Lite™ XLS™+ is the most technically and ergonomically advanced manual multichannel pipette on the market. Equipped with Rainin LiteTouch™ System (LTS™) technology, this multichannel pipette is fast loading, easy on your hand and delivers precise and consistent sample pickup across all channels. Built rugged and with the finest materials available, these multichannels are designed to deliver years of consistent, reliable service. And like all Rainin XLS pipettes, the embedded RFID tag enables advanced calibration tracking.

**METTLER TOLEDO**



# Pipet-Lite XLS+ Multichannel Pipettes

## Performance You Can Feel

### Advanced Functionality

- 8- and 12-channel formats for volumes from 1 to 1200  $\mu\text{L}$
- With LTS for consistent sample pickup on all channels
- Fast easy tip mounting and progressive ejection
- Reliable leak-free seal – every time
- Ergonomic handle design with fingerhook
- Volume adjustment lock prevents drift
- Built-in tip ejector shock absorber reduces thumb stress
- Pipette's RFID tags can be queried for asset management



## Technical Data

### Specifications

Model	Volume ( $\mu\text{L}$ )	Increment ( $\mu\text{L}$ )	Accuracy		Precision	
			%	$\mu\text{L}$ ( $\pm$ )	%	$\mu\text{L}$ ( $\leq$ )
10 $\mu\text{L}$	1.0	0.02	4	0.04	5	0.05
	5.0		1.5	0.075	1.5	0.075
	10.0		1	0.1	0.5	0.05
20 $\mu\text{L}$	2	0.02	7.5	0.15	4	0.08
	10		1.5	0.15	1.5	0.15
	20		1	0.2	0.5	0.1
50 $\mu\text{L}$	5	0.05	3.5	0.18	2	0.1
	25		1.2	0.3	0.5	0.125
	50		0.8	0.4	0.4	0.2
200 $\mu\text{L}$	20	0.2	2.5	0.5	1.2	0.24
	100		0.8	0.8	0.25	0.25
	200		0.8	1.6	0.3	0.6
300 $\mu\text{L}$	30	0.5	2.5	0.75	1.5	0.45
	150		0.8	1.2	0.3	0.45
	300		0.8	2.4	0.25	0.75
1200 $\mu\text{L}$	100	2	3.6	3.6	1	1
	600		0.8	4.8	0.3	1.8
	1200		0.8	9.6	0.2	2.4

### Ordering Information

MT-Order No.	Description	Volume Range
8-Channel		
17013802	L8-10XLS+	0.5–10 $\mu\text{L}$
17013803	L8-20XLS+	2–20 $\mu\text{L}$
17013804	L8-50XLS+	5–50 $\mu\text{L}$
17013805	L8-200XLS+	20–200 $\mu\text{L}$
17013806	L8-300XLS+	20–300 $\mu\text{L}$
17014496	L8-1200XLS+	100–1200 $\mu\text{L}$
12-Channel		
17013807	L12-10XLS+	0.5–10 $\mu\text{L}$
17013808	L12-20XLS+	2–20 $\mu\text{L}$
17013809	L12-50XLS+	5–50 $\mu\text{L}$
17013810	L12-200XLS+	20–200 $\mu\text{L}$
17013811	L12-300XLS+	20–300 $\mu\text{L}$
17014497	L12-1200XLS+	100–1200 $\mu\text{L}$



### Mettler-Toledo Rainin, LLC

Subject to technical changes  
 © 03/2018 Mettler-Toledo Rainin, LLC  
 Printed in Switzerland 17700759 Rev J



RAININ TRADEMARKS: RAININ, Pipet-Lite, XLS, LTS, Pipetting 360° and LiteTouch are trademarks of Mettler-Toledo Rainin, LLC.

[www.mt.com/rainin](http://www.mt.com/rainin)

Für weitere Informationen

